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Orden de Servicio: WO-12073399

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Número de Serie: US924F0850

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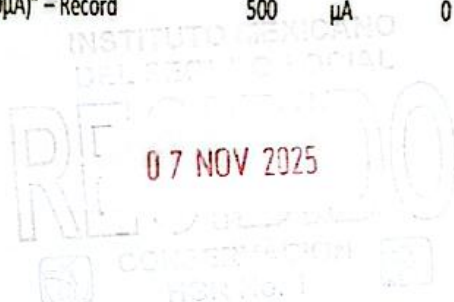
UDI: (01)00884838114180(21)US924F0850

	Min	Max	UoM	Resultado	Conclusión	Comentarios
Planned Maintenance						
System Information						
Software Version				10.0.15		
Hardware Version				B.0		
Inspection & Cleaning						
Query Customer Performance Issues				Pass		
Inspect Exterior Periph. & Power Cords.				Pass		
Inspect & Clean Trackball				Pass		
Inspect & Clean Monitor				Pass		
Inspect/Clean Control Panels & Keyboard				Pass		
Inspected and cleaned the transducer connector assembly (transducer and system side)?				No		
Footswitch, Key- & Remote Control Pad				Not Installed		
Inspect Casters & Brakes Operation				Not Applicable		
Clean/Check System Filters & Fans				Pass		
Check Articulation motion.(monitor, etc)				Pass		
Check System Battery				Pass		
Performance & Peripherals						
Execute Full Set of System Diagnostics				Pass		
Check, Download & Clear Logs				Pass		
Perform Hard-disk Maintenance				Pass		
Backup of System Settings and Configs.				Pass		
Check Operation All User Ctrls. & Lights				Pass		

	Min	Max	UoM	Resultado	Conclusión	Comentarios
Planned Maintenance						
Performance & Peripherals						
Check Operation of Major Applications					Pass	
Test Printers					Not Installed	
Check Operation Dig. Recording Devices					Not Applicable	
Philips Remote Service related Checks						
Is system connected to PRS (Philips Remote Service)?					Yes	
Automatic Export of Log File Configured					Yes	
System Safety General						
Protective Earth resistance (≤ 200 mOhm)	200		mOhm	0	Pass	
Earth Leakage Current. NC. Highest val. (≤ 500 μ A)	500		μ A	0	Pass	
Earth Leakage Current. SF. Highest val. (≤ 1000 μ A)	1000		μ A	0	Pass	
Touch Current 'Normal Condition' (≤ 100 μ A)	100		μ A	0	Pass	
Touch Current – Single fault condition (≤ 500 μ A) – Record higher value between Normal and Reverse	500		μ A	0	Pass	
Transducer Body - BF						
Transducer model					C9-4V	
Transducer serial number					FOLKX2	
Transducer Cables, Connectors & Lens					PASS	
Run Transducer Test					PASS	
Transducer Body - BF						
Transducer model					eL18-4	
Transducer serial number					FOHNDP	
Transducer Cables, Connectors & Lens					PASS	
Run Transducer Test					PASS	
Transducer Body - BF						
Transducer model					C5-1	
Transducer serial number					FOHQ5Y	
Transducer Cables, Connectors & Lens					PASS	
Run Transducer Test					PASS	
Transducer Body - BF						
Transducer model					V9-2	
Transducer serial number					FOGV2W	
Transducer Cables, Connectors & Lens					PASS	
Run Transducer Test					PASS	

PHILIPS

Philips Representative Signature:

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